

CYTOCYPHER



CytoCypher HTS Motorized Microscope System **CYCY100**

User Manual

Copyright 2022 CytoCypher bv

All rights reserved. No parts of this work may be reproduced in any form or by any means - graphic, electronic, or mechanical, including photocopying, recording, taping, or information storage and retrieval systems - without the written permission of the publisher.

Products that are referred to in this document may be either trademarks and/or registered trademarks of the respective owners. The publisher and the author make no claim to these trademarks.

While every precaution has been taken in the preparation of this document, the publisher and the author assume no responsibility for errors or omissions, or for damages resulting from the use of information contained in this document or from the use of programs and source code that may accompany it. In no event shall the publisher and the author be liable for any loss of profit or any other commercial damage caused or alleged to have been caused directly or indirectly by this document.

CytoCypher bv
De Boelelaan 1081
1081 HV
Amsterdam
The Netherlands

<https://www.ionoptix.com/>
info@ionoptix.com
US: +1 617 696 7335
EU: +353 1 685 4800

Document date: November 17, 2022

Research Use Only

This product is intended for research purposes only. It is not certified for clinical applications (including diagnostic purposes). Use of this product in uncertified applications is in violation of FDA regulations.

Table of Contents


Technical Specifications	5
Dimensions and Connections	6
About This Manual	8
Glossary	8
General safety precautions	9
Lid	9
Warnings	9
Product description	9
Installation	10
Disconnecting Cabling	10
Reconnecting cabling:	10
How to use the product	11
Powering-up	11
Temperature control	12
Microscope lamp	12
Switching dichroic cubes inside the microscope	14
Swapping objectives	15
35mm dish holder and perfusion	16
Cleaning and Maintenance	17
Cleaning	17
Maintenance	17
Storage and Transport	17
Storage	17
Transport	17
Faults and Troubleshooting	17
Troubleshooting	17
Repairing	18
Dismantling and Disposing	18




Technical Specifications

Dimensions			
microscope	W x D x H	512 x 274 x 305	mm
controller box	W x D x H	296 x 312 x 133	mm
Weight			
microscope		21	kg
controller box		12	kg
Power			
Voltage		110 – 250 V /	Volt
Frequency		50-60	Hz
Power consumption	Max.	500	VA
Environment			
Operating temperature	Min - Max	15 - 25	°C
Operating Humidity	Min - Max	20 - 80	% RH
Stages			
xy stage	X - Y	18 x 18	mm
z stage	Z	2.5	mm
Resolution xy stage	X - Y	50	nm
Resolution z stage	Z	5	nm
Temperature control			
Setpoint range	Min - Max	25 - 40	°C

Dimensions and Connections

The dimensions of the microscope body (top figures) and controller box (bottom figures). Dimensions are in mm.

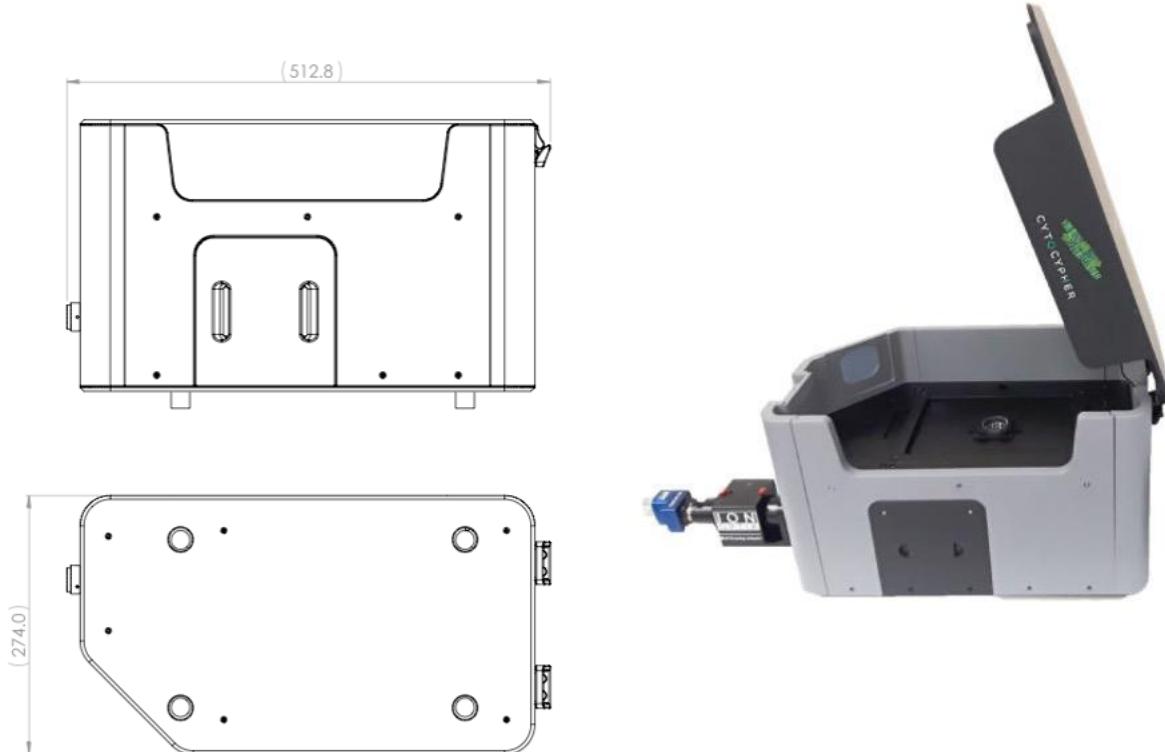

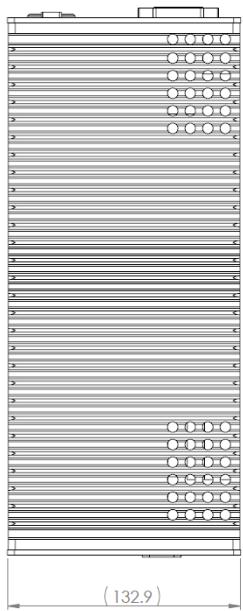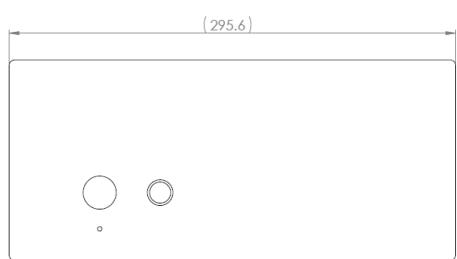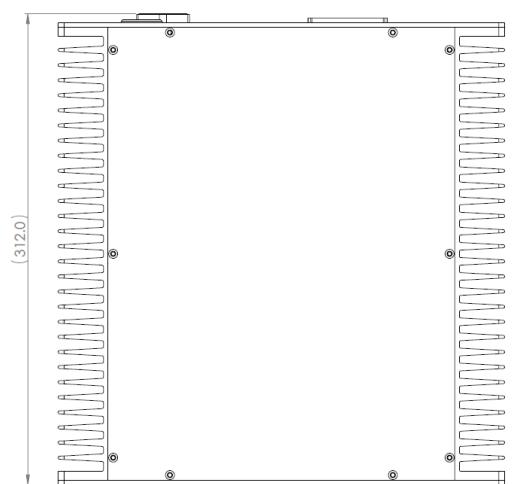

About This Manual

This manual describes how to use the CytoCypher CYCY100 HTS Microscope System.

It is intended as a guide for daily operation and cleaning procedures.

This manual is not intended as a Setup Manual. Installation and setup are included in the delivery of all systems and will be done by IonOptix / CytoCypher engineers.

Glossary

The following terms are used in this manual:

Controller box:

The unit that holds all control electronics for the Microscope. It has the Power button and LED light button on the front panel.

Dish holder:

Drip tray:

Horizontal plate inside the Microscope that collects any spilled liquids and leads them out to the back of the Microscope.

Microscope:

The unit that holds the X-Y-Z stages, light and objective.

Protective hat:

Horizontal plate that mounts on top of the objective inside the Microscope. The Protective hat prevents any spilled liquid from entering the Microscope unit.

Top plate:

Removeable horizontal plate inside the Microscope. This has the mounting position for the dish holder.

General safety precautions

The following warnings and procedures should be adhered to at any time:

- The microscope must be placed on a horizontal surface.
- Use in dry and clean environments only.
- Use only the supplied cables, or cables approved by the manufacturer.

Lid

The microscope has a lid on which the light is mounted.

Before opening the lid, turn off the light with the switch on the front panel of the controller box.

Always open the lid fully so it will stand upright by itself.

Warnings

Always close the lid gently.

Make sure there are no objects or body parts that can get pinched between the lid and the casing.

Product description

The CytoCypher CYCY100 HTS Microscope offers high-speed, high-content mechanical myocyte data acquisition from hundreds of myocytes per hour by using our automated cell finding image analysis method.

The system includes a fast x-y-z position programmable scanning microscope which can mark and find position, size, focus and orientation of myocytes enabling repeated measures of contraction-relaxation function.

Additionally, the enclosed microscope maintains stable temperature while allowing full access for introduction of perfusion systems as well as electrical stimulation leads.

Installation



Attention: The manufacturer will install and setup the device on location.


If the device must be moved to a different location, disconnect all cabling first.

Transport the microscope part and the controller box separately.

Then, reconnect the cabling.

Disconnecting Cabling

- Make sure the power to the device is off, and remove the power plug from the wall-socket
- Disconnect cables between the microscope and its controller box (1 and 2)
- Disconnect the ethernet cable (green) cable to the computer, in the extra card with ethernet slot (1 and 3)
- Make sure the power on the controller box is switched off (switch at front panel).
- Connect the power supply cable.

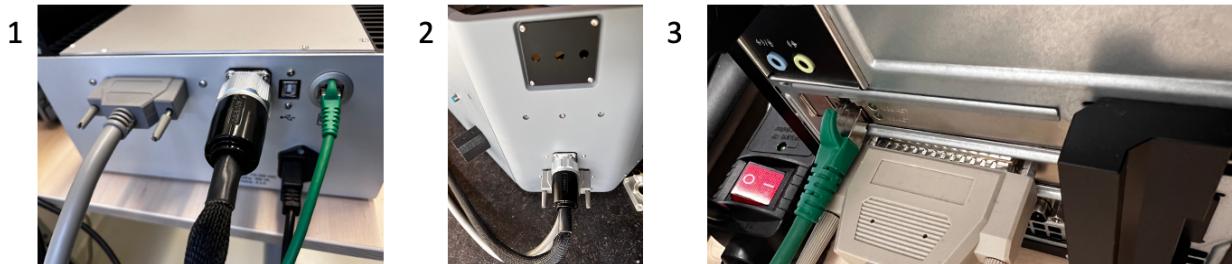

Reconnecting cabling:

- Connect cables between the microscope and its controller box (1 and 2)
- Connect the ethernet cable (green) cable to the computer, in the extra card with ethernet slot (1 and 3)
- Make sure the power on the controller box is switched off (switch at front panel).
- Connect the power supply cable to the controller box and the wall outlet.
- Turn on the system with power switch on the front of the controller box.

How to use the product

Powering-up

To start the system:

- Check that the microscope is empty by opening the lid.
- Loose items can interfere with the movement of the microscope and (if you have one) the Multiwell stage and may damage the system.
- Turn on the system by pressing the power button.
- The indicator light below will be red/orange to indicate the system is initializing. During initialization the objective will move quickly to maximum positions. If the system has a Multiwell stage, during initialization it will move to the upper left corner. When the indicator light turns green, the microscope has initialized successfully and can be used.



Figure 1 Left: system turned off. Right: system turned on and the green light indicates that the system has initialized successfully. Red button on the right is to operate the microscope lamp.



Figure Don't put any dishes in the microscope except in the designated dish holder

Temperature control

The microscope is temperature controlled and by default set to 37°C.

It takes ~30 minutes for the microscope to warm up.

It is advised to power on the microscope at least an hour ahead of the intended start of the experiment.

A small (1-2°C) offset between the temperature of the microscope and the temperature in the dish is expected.

It is advised to occasionally measure the actual temperature in the dish.

The temperature can be controlled in the software (see manual), and at room temperature (20°C) the microscope can be controlled between 25°C and 40°C.

Microscope lamp

The microscope lamp/condenser consists of 2 LED rings, that allows to set up phase contrast (figure 4 and 5).

Switching the Condenser on:

- Switch on the LED ring with the power switch (Figure 3-A) located on the front of the controller box. The switch will light

Switching the Condenser off:

- Switch off the power by pushing the button again (Figure 3-A)



Figure 3 Lamp Power Switch

Change between LED rings

- Turn the LED ring off via the front panel of the Controller Box.
- Open lid and remove the magnetic filter assembly from LED ring condenser (Fig 4)
- Flip switch to the desired position (Figure 5)
- Choose position A / Inner ring when imaging with a 10x Ph1 objective or when imaging with standard objectives.
- Choose position B / Outer ring when imaging with a 20x or 40x Ph2 objective
- Place the filter back onto the condenser.

A **B**

Figure 4 Magnetic filter



Figure 5 LED ring selection. A:
inner ring, B: Outer ring





Switching dichroic cubes inside the microscope

The dichroic cube inside the microscope can be switched to enable fluorescence measurements at different wavelengths.

- Open the port on the side panel of the microscope (Figure 6, left).
- For all the systems delivered after 2020 the port is held shut with magnets, so it can be pulled off; no screws must be removed.
- Reach into the microscope and grab the cube and pull it out (Figure 6, right).
- For all the systems delivered after 2020 the cube is also held with magnets.
- Cubes can be put in only one way, as the magnets will reject every other position. Put the cube back so that it clicks in place.
- Place port cover back on side panel.



Figure 6 Pull of the cover (left) and remove the dichroic (right). Both are held in place with magnets.

Swapping objectives

The microscope can accommodate objectives with an RMS thread and Nikon objectives with a M25 thread. Nikon objectives can be screwed directly into the xyz stage (fig 7). RMS threads need an adapter.

To change an objective, go through the following steps, as depicted in figure 7:

- Turn off the microscope with the switch on the controller box.
- Remove the fixed stage from the microscope by loosening the 4 screws on the corners.
- Loosen the screw under the sticker of the Objective Hat (figure 7.2).
- Remove the Objective Hat (figure 7.3).
- You should now see the objective (figure 7.4).
- Unscrew the objective (figure 7.5).
- Depending on the objective type, you need to use the brass adapter (RMS thread, Olympus in this picture) or, when it has a M25 thread (the Nikon here) without an adapter (figure 7.6).
- After screwing the objective back in, Remove the sticker from the objective Hat and reinstall it (figure 7.7).
- Put a new sticker on the Objective Hat (figure 7.8).

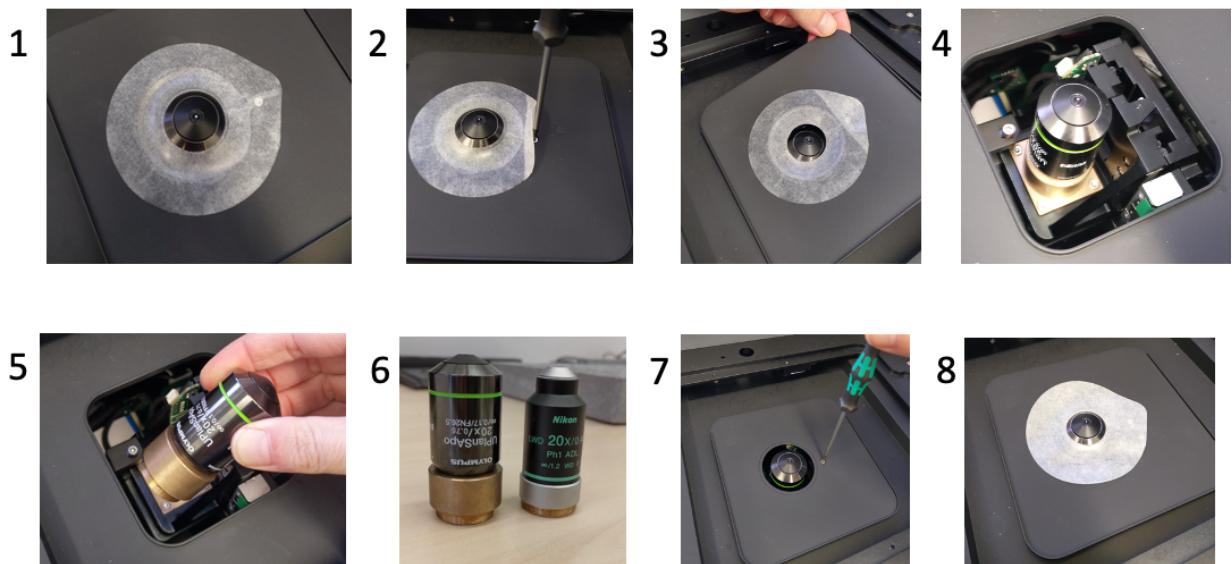


Figure 7 Changing the objective

35mm dish holder and perfusion

The system is equipped with a 35mm dish holder with integrated pacing electrodes and the ability to connect perfusion lines. Figure 8 shows how the dish holder is used.

- 1) The dish holder is placed on the base plate. The holder is dish specific. Please check with the manufacturer for which brands your holder is suitable.
 - 2) The dish holder is placed
 - 3) Fixate the dish holder with the screws provided.
 - 4) place the 35mm dish in the holder. Press down gently so the dish is flush with the bottom
- B1) Place the electrode set on the dish holder. The electrode set is universal for all 35mm dish types.
- 2) The electrode set clicks in place with magnets. It can be placed in 4 different orientations, which one to choose depends on the user. To disconnect, gently pivot electrode set to release the magnets.
 - 3) Guide the wires that connect the pacer with the electrodes set through the ports in the microscope body
 - 4) Connect the wires to the electrode set.
- C1) The electrode set has two holes that fit luer connectors.
- 2) Tubing can be guided through the ports in the microscope body and connected to the electrode set. To perfuse the system, balanced flow is required, i.e. the inflow and outflow should be the same, balanced, in order to prevent flooding the 35mm dish. The CytoCypher MultiPump system can be used for this.

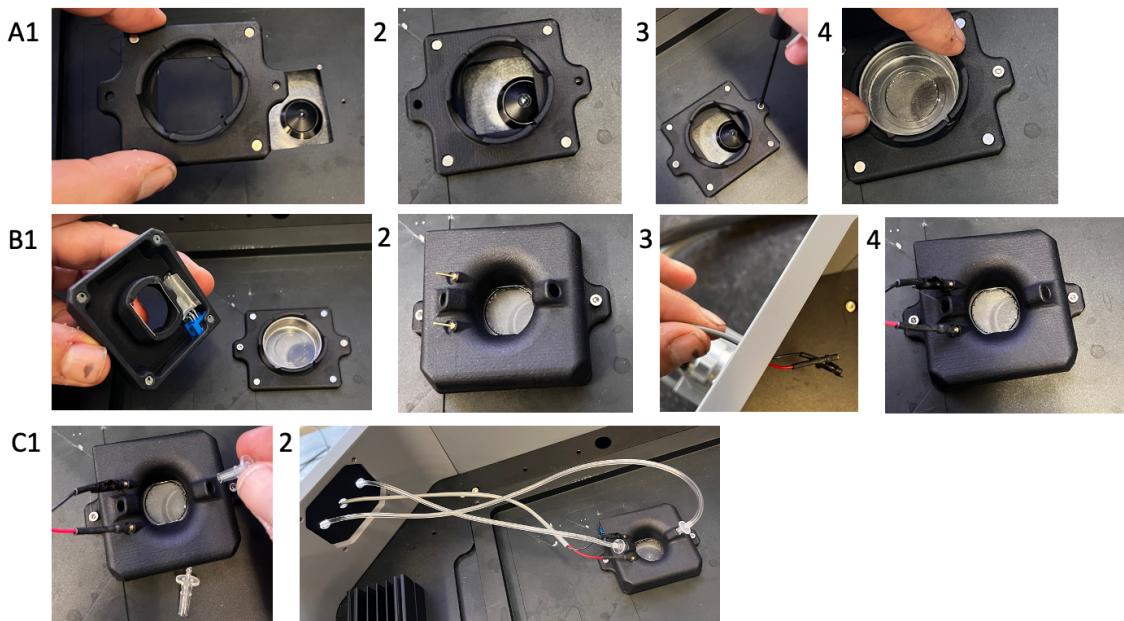


Figure 8 A: placing of 35mm dish holder. B: Placing of electrode set. C: Adding perfusion tubes.

Cleaning and Maintenance

Cleaning

Before cleaning, always turn off the device with the switch on the front panel of the controller box.

Use a soft damp cloth to clean the outside of the microscope and the controller box.

For cleaning the working area in the microscope, 70% alcohol is allowed.

Do not use solvents like acetone or thinner

Make sure the Objective Hat is always present when cleaning, no solution is allowed to drip into the microscope, as it will irreversibly damage the instrument

Maintenance

No maintenance other than cleaning is necessary.

In case of malfunctions contact the manufacturer

Storage and Transport

Storage

When not in use, to protect it from dust, make sure the lid of the microscope is closed.

Store in a dry and clean environment.

Transport

For transport, the original box is preferred.

The high-speed X-Y stage must be locked. Please contact your sales representative for instructions on locking the stage.

Faults and Troubleshooting

Troubleshooting

System does not turn on:

Check if power cord is connected properly to both the controller box and the wall-socket.

Check if power switch on front panel of controller box is in the ON position.

No image / light:

Check that the switch on the front panel of the controller box is in the ON position and the light in the switch lights.

Press the switch to turn on the light.

If the switch is lit, but the light inside the microscope is still not working: check that the switch inside the

light-holder is not in the “OFF” position. (The switch has positions for “Inner ring” / “Outer ring” / “OFF”)

If the light still doesn’t work: check if the cabling between the controller box and the microscope is connected properly and the locking screws are tight.

Repairing

If the device is damaged, stop using it, disconnect all cabling and contact the manufacturer.

The system cannot be repaired by the user, and there are no user-serviceable parts inside.

For repair or inspection, contact the manufacturer.

Dismantling and Disposing

The Multicell microscope is an electronic device and should be disposed of according to local regulations.